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[REF] KED048

ENG ADVANCED WRIST SUPPORT

Please read the following instructions and keep for future reference.

For support, compression and to retain heat around the wrist. Intended uses can include but are not limited to repetitive strain injury (RSI), sprains and instability, carpal tunnel syndrome, tendonitis and arthritis.

The **Kedley Aero-Tech** range is designed for maximum comfort and stability to help recover from injury and prevent further damage. The range has high quality breathable neoprene to retain body heat around the joint and muscles to increase flexibility and stimulate blood flow. The support is lined with a terry-cloth fabric for comfort and to allow airflow. Each item has additional support strategically placed to ensure stability and added strength.

Note: The support should be sized correctly to ensure maximum benefit and to avoid over compression. If the support is being used following an injury, seek medical advice if no improvement in the condition is seen within 72 hours. If the injury is severe or if in any doubt seek medical advice immediately. If used by children, adult supervision is required. Keep out of reach of young children. Always use on clean intact skin. Do not use if liniments, ointments or creams have been applied to the affected area.

Warning: If soreness or a rash occurs discontinue use immediately and seek medical advice. This product contains latex which may cause allergic reactions, do not use if sensitive to latex. *If any serious incidents occur in relation to using this device it should be reported to the manufacturer and the competent authority of the member state in which the user is established.*

Washing & Storage Instructions: Hand wash in lukewarm soapy water, rinse and air dry. Store in a cool, dry place and away from direct heat and sunlight.

AERO-TECH NEOPRENE ADVANCED WRIST SUPPORT

Website: www.kedley.co.uk



Impo International Ltd
Unit 3, 60 Lower Road,
Harrow, Greater London,
HA2 0DH, UK.



Eurolink Europe Compliance Ltd,
Dublin, D02 AY86, Ireland.

REV: KED048IFU-05/21v1

- 1) Check size limit on packaging. Slip thumb through hole ensuring the logo faces outwards.
- 2) Wrap strap down and around wrist. Secure fastener onto the other side of the strap.

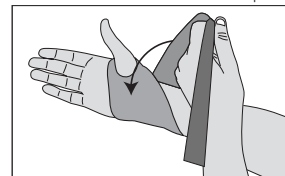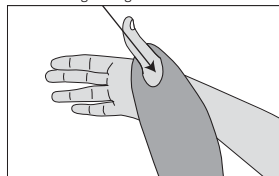- 3) Wrap power band around the wrist for extra compression.
- 4) Check the support fits tightly but does not restrict circulation.

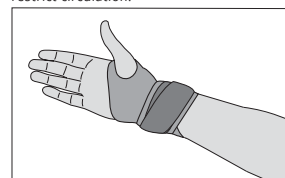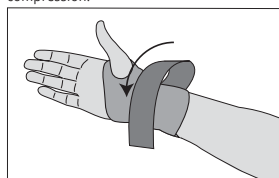KEDLEY

06052021

Video fitting instructions available @ www.kedley.co.uk

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