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ENG HINGED KNEE

Please read the following instructions and keep for future reference.

For support, compression and to retain heat around the knee. Intended uses can include but are not limited to medial/lateral ligament protection, knee instability, knee sprains and strains, hyperextension and meniscus cartilage injuries.

The Kedley Aero-Tech range is designed for maximum comfort and stability to help recover from injury and prevent further damage. The range has high quality breathable neoprene to retain body heat around the joint and muscles to increase flexibility and stimulate blood flow. The support is lined with a terry-cloth fabric for comfort and to allow airflow. Each item has additional support strategically placed to ensure stability and added strength.

Note: The support should be sized correctly to ensure maximum benefit and to avoid over compression. If the support is being used following an injury, seek medical advice if no improvement in the condition is seen within 72 hours. If the injury is severe or if in any doubt seek medical advice immediately. If used by children, adult supervision is required. Keep out of reach of young children. Always use on clean intact skin. Do not use if liniments, ointments or creams have been applied to the affected area.

Warning: If soreness or a rash occurs discontinue use immediately and seek medical advice. If any serious incidents occur in relation to using this device it should be reported to the manufacturer and the competent authority of the member state in which the user is established.

Washing & Storage Instructions: Remove metal hinges then hand wash in lukewarm soapy water, rinse and air dry. Store in a cool, dry place and away from direct heat and sunlight.

AERO-TECH NEOPRENE HINGED KNEE SUPPORT

Website: www.kedley.co.uk


Impo International Ltd
Unit 3, 60 Lower Road,
Harrow, Greater London,
HA2 0DH, UK.


Eurolink Europe Compliance Ltd,
Dublin, D02 AY86, Ireland.

REV: KED025IFU-05/21v1


- 1) Check size limit on packaging. Undo all straps and position patella cut-out over the knee cap.

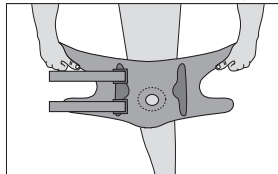
- 2) Wrap the neoprene top and bottom straps around knee and attach.

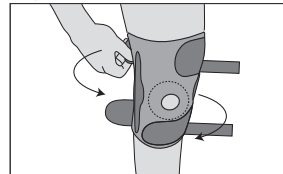- 3) Wrap top and bottom fastening straps around the knee and attach. Re-adjust as required for compression level.

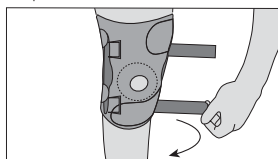- 4) Check the support fits tightly but does not restrict circulation.

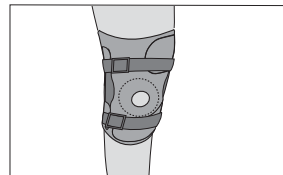KEDLEY

06052021

Video fitting instructions available @ www.kedley.co.uk

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