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[REF] KED045/046

ENG ADVANCED FOAM ARM SLING

Please read the following instructions and keep for future reference.

The **KEDLEY Orthopaedic Soft Neck Collar** has been manufactured using soft foam. The design allows for firm and effective neck support. The adjustable velcro strap enables you to adjust the collar size, whilst ensuring a comfortable fit.

The Kedley Orthopaedic range is made from only the highest quality materials. Each item is manufactured to give you steady support to promote healing and comfort.

Note: The support should be sized correctly to ensure maximum benefit and to avoid over compression. If the support is being used following an injury, seek medical advice if no improvement in the condition is seen within 72 hours. If the injury is severe or if in any doubt seek medical advice immediately. If used by children, adult supervision is required. Keep out of reach of young children. Always use on clean intact skin. Do not use if liniments, ointments or creams have been applied to the affected area.

Warning: If soreness or a rash occurs discontinue use immediately and seek medical advice. This product contains latex which may cause allergic reactions, do not use if sensitive to latex.

If any serious incidents occur in relation to using this device it should be reported to the manufacturer and the competent authority of the member state in which the user is established.

Washing & Storage Instructions: Hand wash in lukewarm soapy water, rinse and air dry. Store in a cool, dry place and away from direct heat and sunlight.

ORTHOPAEDIC SUPPORTS

SOFT NECK COLLAR

Website: www.kedley.co.uk

Impo International Ltd
Unit 3, 60 Lower Road,
Harrow, Greater London,
HA2 0DH, UK.

Eurolink Europe Compliance Ltd,
Dublin, D02 AY86, Ireland.


REV: KED045IFU-05/21v1

KEDLEY


06052021


- 1) Detach the fasteners and open the two ends of the collar.

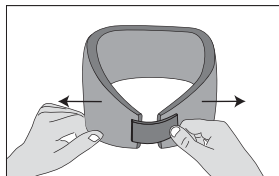
- 2) Look straight ahead and fit the centre of the collar under your chin.

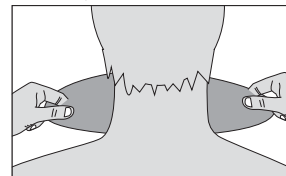- 3) Attach the fasteners at the back of your neck.

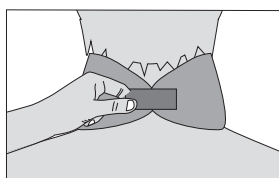- 4) Make sure that you can breathe comfortably. If the collar feels too tight, remove it and put it on again.

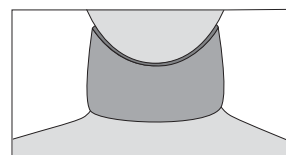Video fitting instructions available @ www.kedley.co.uk

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