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[REF] KED020

ENG ANKLE SUPPORT

Please read the following instructions and keep for future reference.

For support, compression and heat retention around the ankle. Intended uses can include but are not limited to mild ankle sprains and strains, inflammation/swelling, arthritis, soreness, tendonitis, overuse, ankle injuries.

The **PRO-LIGHT** range is made from neoprene rubber and infused with polyester as well as nylon for extra strength, flexibility and comfort, to ensure a secure fit. The neoprene fabric is designed to retain body heat; this warming effect on the muscles and joints increases flexibility and helps prevent injuries. Each product is designed to give maximum support whilst maintaining a comfortable fit.

Note: The support should be sized correctly to ensure maximum benefit and to avoid over compression. If the support is being used following an injury, seek medical advice if no improvement in the condition is seen within 72 hours. If the injury is severe or if in any doubt seek medical advice immediately. If used by children, adult supervision is required. Keep out of reach of young children. Always use on clean intact skin. Do not use if liniments, ointments or creams have been applied to the affected area.

Warning: If soreness or a rash occurs discontinue use immediately and seek medical advice. If any serious incidents occur in relation to using this device it should be reported to the manufacturer and the competent authority of the member state in which the user is established.

Washing & Storage Instructions: Hand wash in lukewarm soapy water, rinse and air dry. Store in a cool, dry place and away from direct heat and sunlight.

PRO-LIGHT NEOPRENE ANKLE SUPPORT

Website: www.kedley.co.uk

Impo International Ltd
Unit 3, 60 Lower Road,
Harrow, Greater London,
HA2 0DH, UK.

Eurolink Europe Compliance Ltd,
Dublin, D02 AY86, Ireland.


REV: KED020IFU-05/21v1


KEDLEY

06052021


- 1) Check size limit on packaging. Slip foot through support ensuring the logo is on the upper strap.

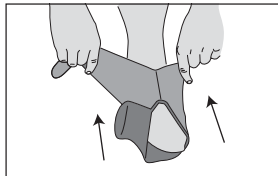
- 2) Secure bottom strap fastener across the bridge of foot.

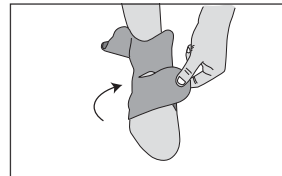- 3) Secure upper strap fastener around ankle.

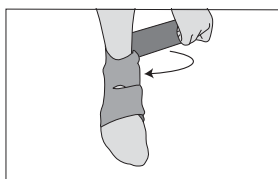- 4) Check the support fits tightly but does not restrict circulation.

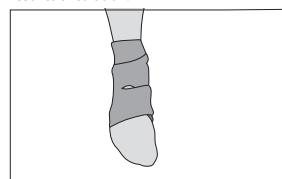Video fitting instructions available @ www.kedley.co.uk

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